

(For scientific research use only, not for clinical diagnosis!)

Canine Stem Cell Factor (SCF) Quantitative Detection Kit (ELISA) Instructions for Use Specification: 96T/48T Catalog Number: SYP-C0024

Purpose: Used to detect canine stem cell factors in serum, plasma, cell culture supernatant and other samples

(SCF) concentration.

Page 1 of 20

Please read the instructions carefully before use. If you have any questions, please contact us via:

Official hotline: 400-999-8863

Technical phone number: 18358180525

Email: UpingBio@163.com

Company website: www.upingbio.com For specific shelf life, please refer to the outer packaging label of the kit. Please use the kit within the shelf life. When contacting us, please provide the product number and production date (see box label) so that we can serve you more efficiently. Page 2 of 20

[Kit performance]

Physical properties: Each liquid component is clear and transparent, with no sediment or floc. Microplate aluminum foil bags should be vacuum packed without damage or leakage.

Calibration curve linearity: The correlation coefficient r value of the calibrator dose-response curve is greater than or equal to 0.9900. Precision: intra-batch variation coefficient CV% is less than 10%; inter-batch variation coefficient CV% is less than 15%. Sensitivity: The lowest detectable dose is less than 7.812 pg/mL.

Recovery rate: The recovery rate is between 85%-115%. Sensitivity: This kit recognizes natural canine stem cell factor (SCF) and has no crossover with structural analogs. Stability: Stored at 2°C-8°C, validity period is 6 months.

Detection range: 31.25 pg/mL - 1000 pg/mL.

Page 3 of 20

[Experimental principle]

This kit uses a double-antibody sandwich enzyme-linked immunosorbent assay (ELISA). In the microwell enzyme plate pre-coated with anti-canine stem cell factor (SCF) antibody (solid-phase antibody), add canine stem cell factor (SCF) calibrator and sample to be tested, then add biotin-labeled antibody, and after incubation and After sufficient washing, HRP-coupled avidin is added. After incubation and sufficient washing, unbound components are removed, and a solid-phase antibody-antigen-biotin-labeled antibody-avidin is formed on the solid surface of the microplate. Enzyme sandwich complexes. Add TMB chromogenic solution to produce a blue product. Under the action of the stop solution, it is finally converted into yellow. The absorbance (OD value) is measured on a microplate reader at a wavelength of 450nm. The absorbance (OD value) is related to the canine stem cell factor (OD value) in the sample to be tested. SCF) concentration is positively correlated. By fitting the calibrator curve, the concentration of canine stem cell factor (SCF) in the sample can be calculated.

Page 4 of 20

[Kit components and storage]

Components		quantity	Main
Calibrator	High Standard	2 vial	Calibrator freeze-
Calibration solution	Reconstitution	2 vial	PBS
Calibrators & Sample Diluents	Standard & Sample Diluent	25mL	PBSTN
coated microplate	Microelisa Stripplate	96T/48T	Pre-coated solid phase
biotin antibody	Bio-Antibody	10mL	biotin antibody
HRP labeled avidin	HRP- Conjugate	10mL	HRP labeled avidin
TMB chromogenic	TMB	10mL	TMB
stop solution	Stop Solution	6mL	acidic solution
20×concentrated	20X Wash Solution	25mL	0.05%Tween20
manual	manual	1 serving	
Ziplock bag	Ziplock bag	1	
Self-adhesive	Self-adhesive	4 pieces	

Note: 1. Before use, please check whether the label and quantity of the reagents in the kit are consistent with the table.

2. The test kit should be stored at 2-8°C. Expired test kits must not be used.

3. If the coated microplate is not used up in a single time, remember to seal it and store it at 2-8°C.

4. The reconstituted calibrator can only be used on the same day.

5. If the components of the kit need to be used again, please ensure that they have not been contaminated since the last use.

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 5 of 20

UpingBio technology Co.,Ltd

Prepare your own test equipment required for the

test (not provided, but can assist in purchasing) 1.

Standard specification microplate reader.

2. Automatic plate washing machine.

3. Oscillator.

4. A series of adjustable pipettes and tips. When testing a large number of samples at one time, it is best to use a multi-channel pipette.

[Kit limitations]

1. For scientific research use only and not

for clinical diagnosis.

2. Use within the validity period marked on the kit. Expired products must not be used.

3. Do not mix with kits or components from other manufacturers.

4. Use the sample diluent provided with the kit.

5. If the sample value is higher than the highest calibrator concentration value,

please dilute the sample appropriately and then re-measure. 6. Human anti-mouse

and other heterophilic antibodies present in the sample to be tested will interfere with the test results. Please eliminate this factor before testing.

7. The test results obtained by other methods are not directly comparable to the test results of this kit.

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 6 of 20

(Precautions)

1. This kit is for in vitro research only and not for clinical diagnosis.

Please wear a lab coat and latex gloves for protection during the test.
Especially when testing blood or other body fluid samples, please follow the national biological laboratory safety protection regulations.

3. Incubate strictly according to the specified time and temperature to ensure accurate results. All reagents must reach room temperature 20-25°C before use. Store reagents refrigerated immediately after use.

4. Incorrect plate washing may lead to inaccurate results. Be sure to absorb as much liquid as possible from the wells before adding substrate. Do not allow the microwells to dry out during incubation.

5. Eliminate residual liquid and fingerprints on the bottom of the plate, otherwise it will affect the OD value.

6. The substrate chromogenic solution should be colorless or very light in color.

7. Avoid cross-contamination of reagents and specimens to avoid erroneous results.

8. Avoid direct exposure to strong light during storage and incubation.

9. Open the sealed bag after balancing to room temperature to prevent water droplets from condensing on the cold slats.

10. Any reaction reagents must not come into contact with bleaching

solvents or strong gases emitted by bleaching solvents. Any bleaching

ingredients will destroy the biological activity of the reagents in the kit.

11. The microplate reader used for detection needs to be equipped with a filter capable of detecting a wavelength of 450 ± 10 nm.

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 7 of 20

Optical density ranges from 0-3.5. It is recommended to preheat 15 minutes in advance before use.

12. Do not mix or replace the reagents in this kit with reagents from other batch numbers or other sources. 13. The EP tubes and suction tips used in the test are single-use and are strictly prohibited from mixing.

14. Do not use expired reagents.

[Preparation and preservation of samples]

The following lists only general guidelines for sample collection and preservation. During the collection and storage of all samples, sodium azide shall not be used as a preservative. If the sample is not analyzed immediately, it should be aliquoted and stored frozen, and repeated freezing and thawing should be avoided.

Cell culture supernatant: Centrifuge to remove precipitate, analyze immediately or aliquot and store frozen at -20°C.

Serum: Collect blood in a clean test tube, coagulate at room temperature for 30 minutes, centrifuge at 2000×g for 20 minutes, and collect serum. Analyze immediately or aliquot and store frozen at -20°C. Plasma: Anticoagulate with heparin, citrate or EDTA, and centrifuge at 2000 × g for 20 minutes at 2-8°C within 30 minutes of blood drawing. To eliminate the influence of platelets, it is recommended to further centrifuge at 10,000 × g for 10 minutes at 2-8°C. Analyze immediately or aliquot and store frozen at -20°C.

Cell lysis medium: For adherent cells, remove the culture medium and replace with PBS, physiological saline or blood-free

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 8 of 20

UpingBio technology Co.,Ltd

Wash with culture medium. Add an appropriate amount of lysis solution and pipet several times to ensure full contact between the lysis solution and cells. Typically after 10 seconds, cells are lysed. For suspended cells, collect the cells by centrifugation and wash them once with PBS, physiological saline or serum-free culture medium. Add an appropriate amount of lysis solution, pipet with a pipette to disperse the cells, and flick with your fingers to fully lyse the cells. After full lysis, centrifuge at 10000-14000×g for 3-5 minutes and take the supernatant. Analyze immediately or aliquot and store frozen at -20°C.

Urine: Collect in sterile tubes and centrifuge at 2000×g for 20 minutes. Carefully collect the supernatant. If a precipitate forms, centrifuge again.

[Reagent preparation]

1. Before use, all components must be rewarmed for at least 120 minutes to ensure full rewarming to room temperature.

2. Concentrated washing liquid: The concentrated washing liquid taken out from the refrigerator will produce crystals. This is a normal phenomenon. Heating in a water bath will completely dissolve the crystals. Concentrated detergent and distilled water, dilute 1:20, that is, 1 part of concentrated detergent, add 19 parts of distilled water.

Website: www.upingbio.com

Official hotline: 400-999-8863


Supervision phone number:

Page 9 of 20

[Calibration dilution method]

校准品重溶步骤:用校准品复溶液重溶校准品,将一瓶校准品复溶液 中的液体,全部加入到一瓶的校准品冻干粉中,轻柔地涡旋震荡,确保充 分混匀,重溶后校准品的母液浓度为2000 pg/mL ,稀释前充分混匀。

校准品母液稀释步骤:稀释前校准品工作液静置 1-2 分钟,用校准品 & 样本通用稀释液将校准品母液进行倍比稀释,倍比稀释方法:取7支 EP 管,每管中加入 500 µL 校准品&样品稀释液,从 2000 pg/mL 的校准 品 母液中吸取 500 µL 到第一支 EP 管中混匀配成 1000 pg/mL 的校准品工 作液,按此步骤往后依次吸取混匀。如下页图示。

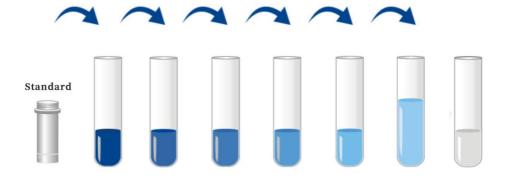

2000 pg/mL

校准品建议稀释浓度:建议配制成以下浓度:1000、500、250、125、

网址: www.upingbio.com

官方热线: 400-999-8863

监督电话: 18158115141

第10页共20页

UpingBio technology Co.,Ltd

62.5、31.25、0 pg/mL,并作为拟合标曲的校准品浓度值。

提示:最后一管直接加入校准品&样本通用稀释液作为0值,不需要 再从倒数第二管中吸取液体。倍比稀释的校准品工作液需要现配现用。

【操作程序】

推荐样本稀释方案:建议老师先做预实验摸索样本最佳稀释倍数,然 后再做正式实验。

所有试剂和组分都先恢复到室温,校准品、质控品和样品,建议做复 孔。

1、按前面说明书描述的方法, 配制好试剂盒各种组分的工作液。

2、从铝箔袋中取出所需板条,剩余的板条用自封袋密封放回冰箱。

3、设置校准品孔、样本稀释液孔、空白孔和样本孔,校准品孔各加 不同浓度的校准品 50μL,样本稀释液孔加样本稀释液 50μL,空白孔不 加,样本孔加待测样本 50μL。除空白孔外,每孔加入生物素抗体 100uL,用封板膜盖住反应板,37℃水浴锅或恒温箱避光温育 60min。

4、揭开封板膜,弃去液体,吸水纸上拍干,每孔加满洗涤液,静置 ^{网址:} www.upingbio.com ^{官方热线:} 400-999-8863 监督电话: 18158115141

UpingBio technology Co.,Ltd

1 minute, shake off the washing liquid, pat dry on absorbent paper, repeat this 5 times. If using an automatic plate washer, please wash the plate according to the operating procedures of the plate washer. Adding a soaking program for 30 seconds can improve the detection accuracy. After washing the plate and before adding substrate, pat the reaction plate dry on clean, lint-free paper.

5. Except for the blank wells, add 100uL of HRP-labeled avidin to each well, cover the reaction plate with a sealing film, and incubate in a 37°C water bath or incubator in the dark for 20 minutes.

6. Repeat step 4.


7. Add 100μL of TMB chromogenic solution to all wells. Cover the reaction plate with sealing film and incubate in a 37°C water bath or incubator in the dark for 15 minutes.

8. Add 50 μ L of stop solution to all wells, and read the absorbance (OD value) of each well on a 450 nm wavelength microplate reader.

Page 12 of 20

[Result calculation]


1. Use the concentration of the calibrator as the abscissa and the corresponding absorbance (OD value) as the ordinate. Use computer software and four-parameter Logistic curve fitting (4-pl) to create a standard curve equation. Through the absorbance (OD value) of the sample value), use the equation to calculate the concentration value of the sample. [Calculation using ELISA Calc software] 2. If the sample is diluted, the concentration value measured by the above method must be multiplied by the dilution factor to obtain the final concentration of the sample.

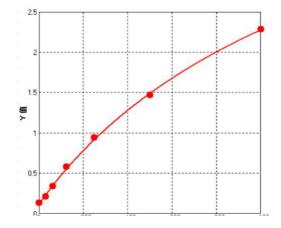

(Schematic diagram, for reference only)

Page 13 of 20


(Operation Summary **)**

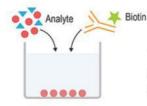

1、反应板孔中加入50uL校准品工作液 或样本后,立即每孔加入100uL生物素化 抗体工作液,37℃孵育60分钟。



2、弃掉板内液体,洗板5次。

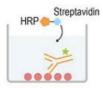

3、每孔加入100uL HRP酶结合物工作液 37℃孵育20分钟,弃掉板内液体,洗板 5次。



4、每孔加入100uL TMB显色液, 37℃孵 育 15分钟。



5、每孔加入 50uL 终止液。



6、立即在 450nm波长下读数,处理数据。

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 14 of 20

(problem analysis **)**

Problem Description	Possible Causes	Corresponding	
	Incorrect liquid	Check pipettes and tips	
Negative and positive	Equilibration time is too	Ensure sufficient	
control results are unstable	Incomplete washing	Ensure the washing time and number of washes	
	Incubation time too short	Ensure adequate	
	The experimental	Use recommended	
Very weak or colorless	Insufficient reagent	Check the liquid	
	Incorrect dilution	aspiration and addition	
	Enzyme label	Mix enzyme conjugate and substrate and check	
Reading value is low	Microplate reader	Check the wavelength and filter settings on the	
	settings are incorrect	Turn on the microplate	
Large coefficient of	Adding fluid incorrectly	Check the filling	
	The working	Use the recommended	
High background value	Incomplete washing of enzyme plate	Ensure that each step of cleaning is complete; if using an automatic plate washer, please check	
	The lotion is	Prepare new lotion	
Low sensitivity	Improper storage of	Store relevant reagents according to instructions	
	Not terminated before	Stop solution should be added to each well	

UpingBio technology Co.,Ltd

If the experimental results are not good, please take pictures of the color development results in time, save the experimental data, keep the strips used and unused reagents, and then contact our company's technical support to solve the problem for you.

[statement]

1. Due to existing conditions and scientific and technological level, it is not possible to conduct comprehensive identification and analysis of all raw materials. This product may have certain quality and technical risks.

 This kit removes/reduces some endogenous interfering factors in biological samples during the development process. Not all possible influencing factors have been removed.

3. The final experimental results are closely related to factors such as the effectiveness of the reagents, the relevant operations of the experimenter, and the experimental environment at the time. Our company is only responsible for the kit itself and is not responsible for the sample consumption caused by the

use of the kit. Please use The user should fully consider the possible usage of the sample and reserve sufficient samples before use.

4. In order to achieve good experimental results, please only use the reagents provided in our company's kits, do not mix products from other manufacturers, and operate in strict accordance with the instructions.

5. Due to incorrect reagent preparation and microplate reader parameter settings during the operation, abnormal results may occur. Please read the instructions carefully and adjust the instrument before the experiment.

Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

Page 16 of 20

UpingBio technology Co.,Ltd

6. Even if operated by the same personnel, different results may be obtained in two independent experiments. In order to ensure the reproducibility of the results, it is necessary to control every step of the experimental process.

7. The kits will undergo strict quality inspection before shipment. However, due to factors such as transportation conditions, differences in experimental equipment, etc., user test results may be inconsistent with factory data.

8. This kit has not been compared with similar kits from other manufacturers or products using different methods to detect the same target, so inconsistent test results cannot be ruled out.

9. The kit is for research use only. If it is used for clinical diagnosis or any other purpose, our company will not be responsible for any problems arising therefrom, nor will we assume any legal liability.

Page 17 of 20

UpingBio technology Co.,Ltd



Website: www.upingbio.com

Official hotline: 400-999-8863

Supervision phone number:

[Experimental experience]

[Experimental experience]